





## Gravity Bag: Tube Feeding Instructions

1. Wash your hands.
2. Gather your supplies in a clean and comfortable place. You will need:
  - Formula
  - Gravity feeding bag
  - IV pole or hook for hanging the feeding bag
  - Cup of water
  - 60 mL catheter tip syringe
3. Close the roller clamp on the feeding bag.

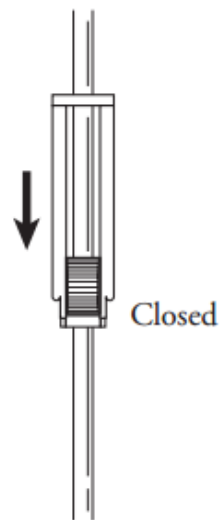


4. Pour the prescribed amount of formula into the feeding bag.

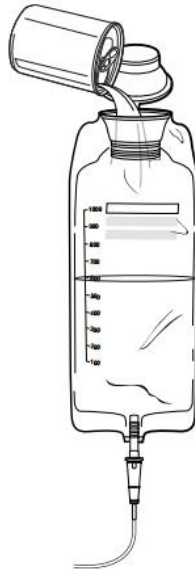




5. Hang the feeding bag above you on either an IV pole or hook. It should be close to a comfortable chair or your bed.

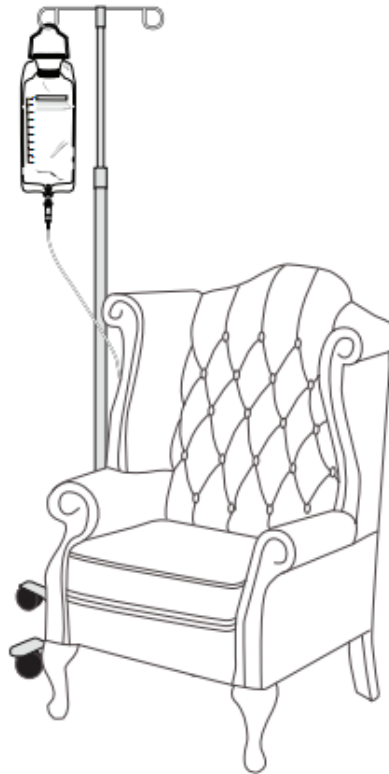

6. Hold the end of the feeding bag tube over an empty cup or the empty formula container. Remove the cap at the end of the tube. Slowly open the roller clamp on the feeding bag tube. Allow the formula to run through the tubing until it reaches the bottom of the tubing, then close the clamp. This gets the air out of the tubing.

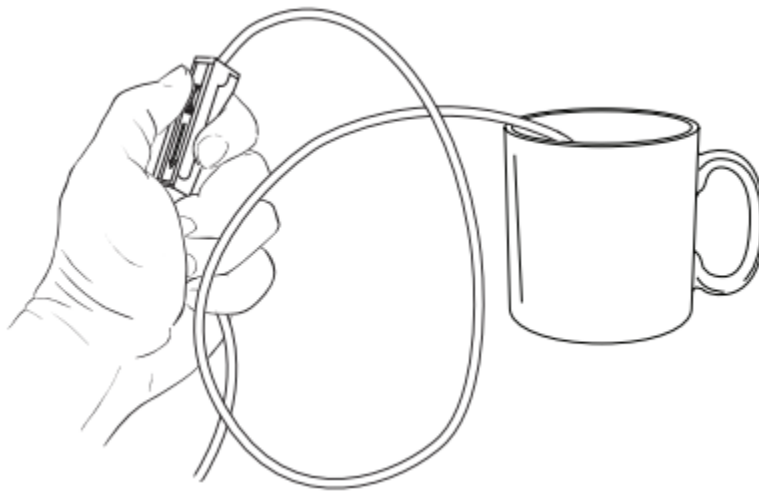




7. Fill the syringe with 30-60cc of water. Place the syringe into your feeding tube and flush your feeding tube.




8. Remove syringe and connect the end of the feeding bag tubing to your feeding tube.

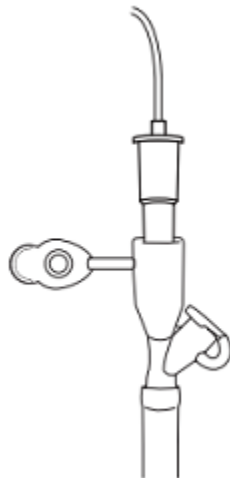



9. Slowly open the roller clamp on the feeding bag tubing. Allow the formula to drip in over 20 to 45 minutes. If the formula is flowing too fast or too slowly, adjust the roller clamp as needed. You can also lower or raise the bag to change the flow of the formula.

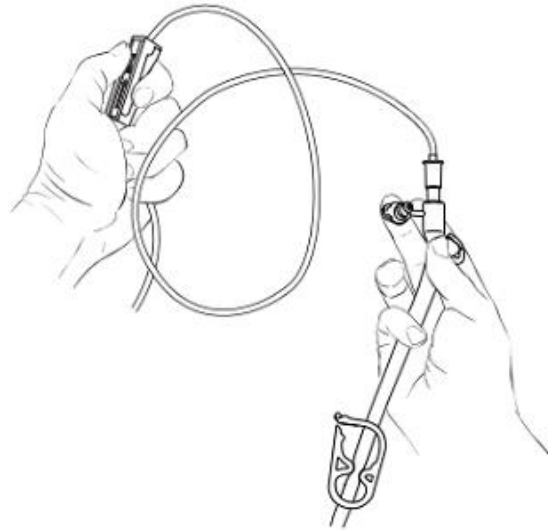

10. When the feeding bag is empty, close the roller clamp on the feeding bag tubing. Disconnect the feeding bag and set it aside.

11. Fill the syringe with 30-60cc of water. Place the syringe into your feeding tube and flush your feeding tube. Close cap of feeding tube when done.



12. Rinse the feeding bag and syringe with warm water after each feeding. Let your supplies air dry. Use a new feeding bag every 24 hours to prevent infection.